

The Integris[®] B Series Headwall Systems



The Integris B® Series Headwall Systems



The Integris® B Series Headwall Systems are patient care in action. These highly functional headwall systems adapt to various patient environments. With appealing lines, the attractive headwall systems help create a pleasant atmosphere in patient rooms and, as the acuity level of patients under care changes, the characteristics of the Integris® B Series Headwall Systems respond.

The streamlined Integris® B Series Headwall Systems are designed to integrate with beds, stretchers, communication systems and other Hill-Rom® equipment. Available in a wide array of High Pressure Laminate (HPL)


surfaces to complement room design and color scheme; its sculpted appearance not only softens the room, but minimizes dust collection and facilitates cleaning.

As with all patient-care environment products at Hill-Rom, our specialists are there to carefully analyze the clinical needs and the unique culture of your unit. Let our team of skilled professionals in architecture, patient beds and surfaces, maternal-infant care, communications, health care furnishings, medical gases and finance help unfold the possibilities. The result... facilities and equipment that blend leading-edge technology within an environment of care.

The Integris® 2002B Headwall System

The Ultimate Headwall System for High Acuity Areas.

- 1 Moveable Gas Outlets allow optimal positioning and serviceability
- 2 Standard 14-Port Manifold located behind gas panel (maximum of 4-Oxygen Outlets, 4-Medical Air Outlets and 6-Vacuum Outlets)
- 3 Optional integrated Bed Locator protects walls and equipment
- 4 Seamless horizontal equipment rail
- 5 Nurse Call and other low voltage provisions

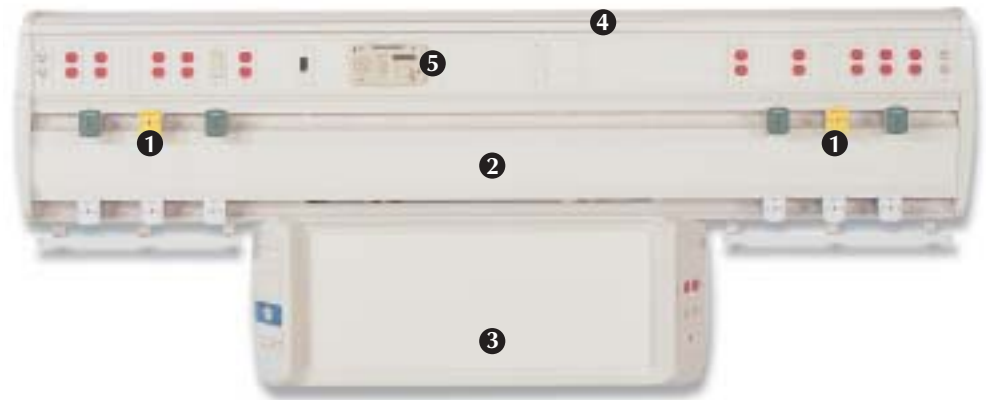

The Integris® 2002B Headwall System organizes medical services conveniently within reach of caregivers. Its exclusive moveable gas outlets allow caregivers to quickly configure the headwall system to address patient needs and staff preference. And, the availability of multiple rails allows for versatile positioning of electrical receptacles, low voltage devices and accessories.

The Integris® B Series Headwall Systems are simpler to install than in-wall services, and the chased units can be mounted on existing walls in any care area to minimize renovation time and construction costs. Additional gas outlets may also be installed on the Integris® 2002B Headwall System after initial installation.

The Integris® 2001B Headwall System

Economy, Flexibility, Dependability and Productivity.

- 1 Fixed Gas Outlets can be located along the rail
- 2 Seamless horizontal equipment rail
- 3 Optional Accessory Bar provides second level of equipment management
- 4 Shown as a private unit. Also available with or without a chase and in semi-private configurations



A more affordable solution that still offers many of the features of the Integris® 2002B is the Integris® 2001B Headwall System. For acuity areas

that require centralization and organization of services, this system can still be customized to meet immediate and future care needs.

The Integris B® Series Dimensions

	INTEGRIS® 2001B	INTEGRIS® 2002B	INTEGRIS® 2003B
Overall			
Length (1 rail)	6'6" - 9'6" (78 - 289.6 cm)	6'6" - 9'6" (78 - 289.6 cm)	6'6" - 9'6" (78 - 289.6 cm)
Rail Height	14.767" (37.5 cm)	24.688" (62.7 cm)	9.451" (24 cm)
Rail Depth	4.40" (11.2 cm)	4.40" (11.2 cm)	10.75" (27.3 cm)
Device Centerline from Floor	42.75" (108.6 cm)	42.75" (108.6 cm)	69.5" (176.5 cm)
Vertical Chase			
Width	13.25" (33.7 cm)	13.25" (33.7 cm)	13.25" (33.7 cm)
Depth	5.5" (14 cm)	5.5" (14 cm)	5.5" (14 cm)
Height	varies with ceiling height	varies with ceiling height	varies with ceiling height

DIRECT CUSTOMER SERVICE HOTLINE:
USA 800-445-3730
CANADA 800-267-2337

Hill-Rom®

A HILLENBRAND INDUSTRY

Batesville, IN 47006 USA • 812-934-7777 • Fax 812-934-7038

Canada • 905-206-1355 • Fax 905-206-0561

International • 812-934-8173 • FAX 812-934-7191

www.hill-rom.com

CTG223ra



Nursemate cabinet (as shown on cover)
is no longer available.

Hill-Rom reserves the right to make changes without notice in design, specifications, and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.
©2002 Hill-Rom Services, Inc.
ALL RIGHTS RESERVED.